

Bed Reference Guide

Glossary

BED WIDTH AND LENGTH

Standard senior living beds are 36" wide. A wider 39", 42", or 48" bed provides more space to sleep and reposition, which increases comfort and could reduce falls. They may feature:

- A specific, fixed width and/or length
- · Built-in adjustments allowing for multiple width options
- Detachable accessories that allow the bed to expand to other widths
- Adjustable lengths

HEIGHT/TRAVEL RANGE

Indicates both the low height and high height of the bed

DECK STYLE

Good decks provide a firm sleeping surface that won't change over time. The most common styles:

Metal slats provide a strong, supportive frame and are easy to clean.

Large-diameter wire welded to a frame provides high mattress ventilation and eliminates sharp edges.

Usually consists of a thick wire grid attached to a metallic frame by springs.

Pan

Offers the best orthopedic support, has no sharp edges, is easy to clean, and comes in two styles:

- · Perforated, which allows for high mattress ventilation
- Composite, which allows for ease in removal and cleaning

LOCKING MECHANISMS

Step Lock

An easy-to-use safety lock that can be initiated by the foot to provide immediate immobilization of the bed.

Floor Lock


Typically extended to place the rubber foot of the lock onto the floor surface to create drag or friction to immobilize the bed. The placement and number of locks is determined by the caster type.

Standard Caster Locks

These locks are individually applied to each caster.

AUTOMATIC CONTOUR

Activated by simultaneously pressing the head rise/lower and foot rise/lower buttons on the bed control to initiate a head and foot cradle effect. Some bed controls have a default automatic contour feature triggered by pressing one button.

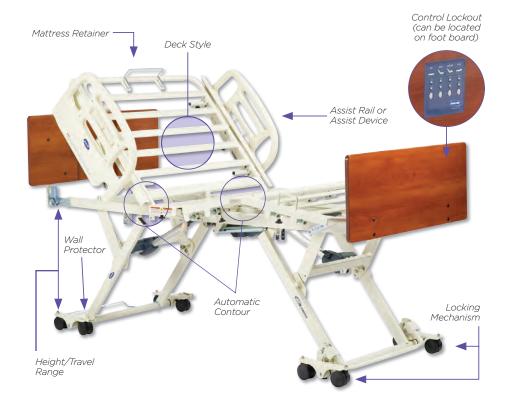

CONTROL LOCKOUT

A resident safety feature that allows the caregiver to lock the hand pendant functionality.

MOBILITY

Enables bed mobility via the bed's wheel locking mechanism, depending on the high/low position of the bed. For example, some beds are only mobile when the bed is in its lowest position, while others are mobile at any height.

TRENDELENBURG POSITION

The Trendelenburg position is a patient position that is used for certain kinds of medical procedures. When a patient is placed in this position, he or she is lying face up, and the surface on which the patient lies is angled so that the feet are above the head. In the reverse position, the head lies above the feet.

See our full-line Catalog V64, pages 1714-1721, and our Healthcare Catalog V9, pages 122-129. for our complete selection of beds.



© 2014 HDS IP Holding, LLC. All Rights Reserved. For a full description of our terms and conditions, see the back of Catalog V64, page 362 in Healthcare Catalog V9, or visit hdsupplysolutions.com/terms.



SUPPLY Bed Reference Guide

Special Orders

BRAND	PRODUCT NAME	WEIGHT CAPACITY	DECK LENGTH	DECK WIDTH	ADJUSTABLE LENGTH	TRAVEL RANGE	DECK Style	LOCKING MECHANISM	LOCKS AT Any Height	INCLUDES H/F BOARDS	INCLUDES RAILS	AUTOMATIC CONTOUR	# OF MOTORS	CONTROL LOCKOUT AVAILABLE	MOBILITY AT Any Height	WARRANTY	PART #
INVACARE Yes, you can.º	CS9 FX600	600 Lbs	80"	36, 39 Or 42"	No	9 To 30"	Slat	Step Lock	Yes	Yes	No	Yes	3	Yes	Yes	5-Year Mechanical 3-Year Warranty 2-Year Electronics	677420, 677421, 677422, 677423, 677425, 677426, 677427, 677428
	CS7	500 Lbs	76 To 80"	36"	Yes	7 To 30"	Slat	Floor Lock System	Yes	No	No	Yes	3	Yes	Yes	10-Year Frame 5-Year Mechanical 3-Year Warranty 2-Year Electronics	676942
	CS5	500 Lbs	76 To 80"	36"	Yes	9 To 30"	Slat	Step Lock	Yes	No	No	Yes	3	Yes	Yes	10-Year Frame 5-Year Mechanical 2-Year Electronics	676939
	CS3	450 Lbs	76 To 80"	36"	Yes	8-1/2 To 27"	Slat	Maxi-Lock Floor Lock System & Standard Caster Locks	No	No	No	Yes	3	Yes	No, Only In Lowest Position	10-Year Frame 5-Year Mechanical 2-Year Electronics	676936
	SC900DLX	450 Lbs	76 To 80"	36"	Yes	9-7/8 To 20-1/4"	Grid	Standard Caster Locks	Yes	Yes	Yes	No	3	No	No, Only In Lowest Position	10-Year Frame And Welds 2-Year Mechanical 1-Year Electronics	673507, 673508
	820-3MDXL	450 Lbs	76 To 80"	36"	Yes	13-1/2 To 27"	Grid	Standard Caster Locks	Yes	Yes	Yes	No	3	No	No, Only In Lowest Position	10-Year Frame And Welds 2-Year Mechanical 1-Year Electronics	673502
NO DA	Elite Riser	600 Lbs	76, 80 And 84"	35"	No	7-1/2 To 28"	Slat	Step Lock	Yes	No	No	Yes	3	Yes	Yes	Lifetime Welds 4-Year Mechanical 4-Year Electronics	400561, 400562, 400563
♦Joerns	UltraCare XT	500 Lbs	76 To 80"	36"	Yes	7 To 30"	Grid	Step Lock	Yes	No	No	Yes	3	Yes	Yes	Lifetime Welds 15-Year Frame 5-Year Parts And Components	679389
	UltraCare XT	500 Lbs	76 To 80"	36"	Yes	7-15/16 To 30"	Grid	Step Lock	No	No	No	Yes	3	Yes	No, Only In Lowest Position	Lifetime Welds 15-Year Frame 5-Year Parts And Components	679388
	EasyCare	500 Lbs	76 To 80"	36"	Yes	7 To 30"	Grid	2-Pedal Locking System	Yes	No	No	No	3	Yes	Yes	Lifetime Welds 15-Year Frame 4-Year Parts And Components	679410
	EasyCare	500 Lbs	76 To 80"	36"	Yes	7 To 30"	Grid	2-Pedal Locking System	Yes	No	No	No	3	No	Yes	Lifetime Welds 15-Year Frame 4-Year Parts And Components	679409
drive	Long Term Care Bed	450 Lbs	80"	36"	No	15 To 26"	Spring Deck	Standard Caster Locks	Yes	Yes	Yes	Yes	3	No	Yes	Limited Lifetime Welds And Frames 2-Year Electronics 1-Year Parts And Components	668059
GENDRON Bariatric Rehab Leadership	BariLow EC	850 Lbs	84"	36, 42, Or 48"	No	9-3/4 To 26"	Pan	Standard Caster Locks	Yes	Yes	Yes	Yes	3	No	Yes	5-Year Frame 1-Year Parts And Components	679296
	Maxi Rest	750 Lbs	80"	48"	No	15-1/2 To 29-12"	Pan	Standard Caster Locks	Yes	Yes	Yes	Yes	3	No	Yes	5-Year Frame 1-Year Parts And Components	676428